



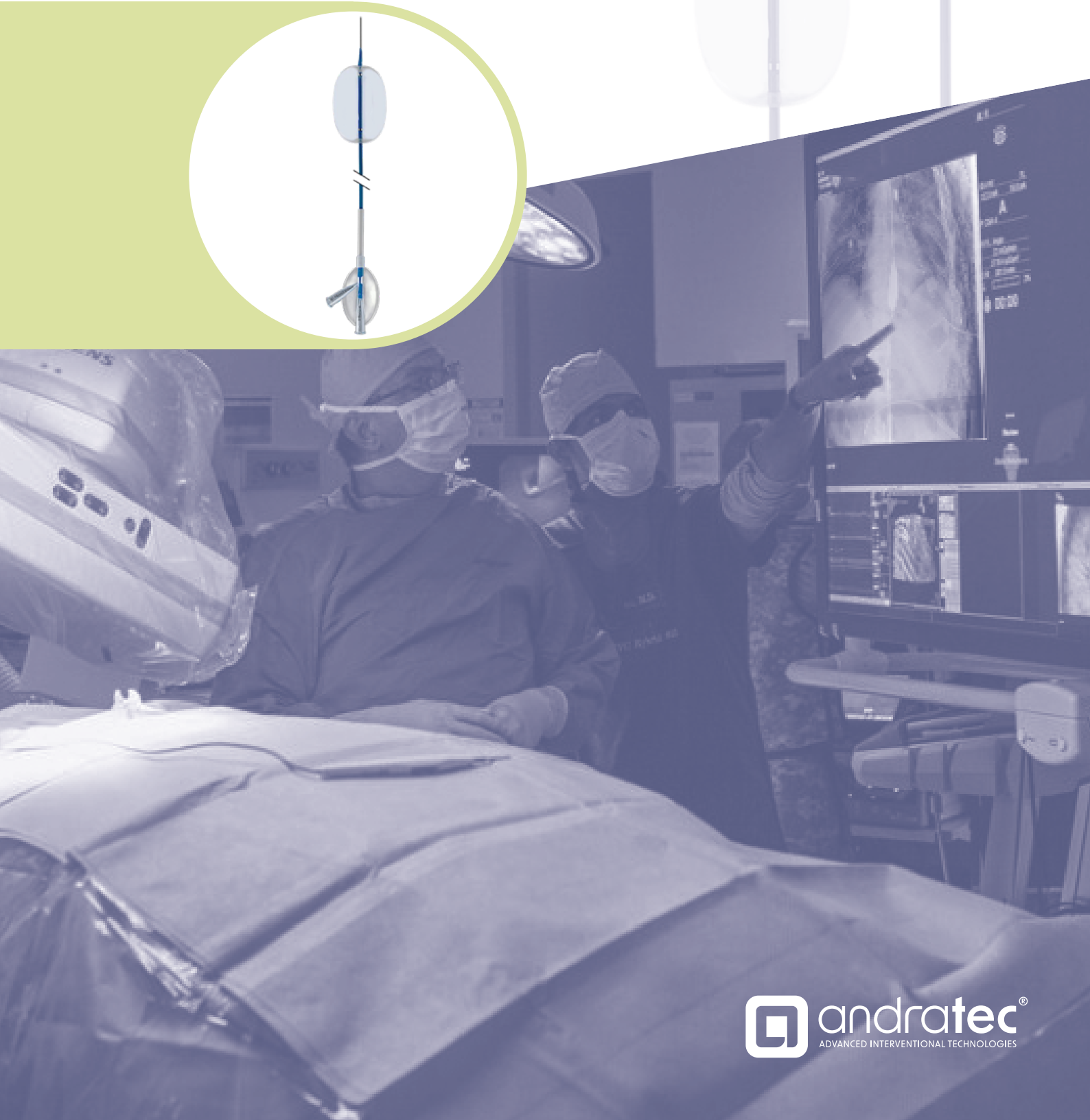

altosa®
SFT | BALLOON
CATHETER

compliance in tacking peripheral stent grafts

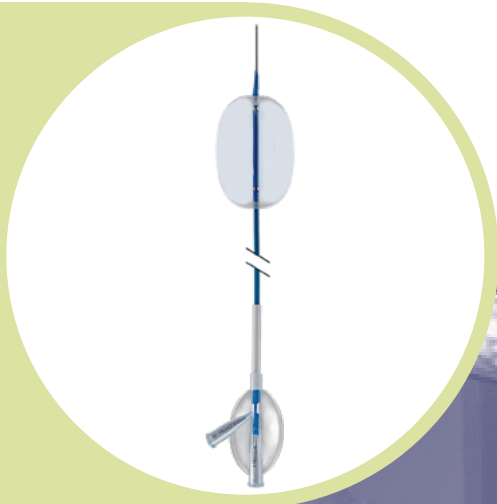

- Especially designed for large lumen dilatation and/or tacking of peripheral endoprosthesis.
- Enhanced inflation/deflation times for safe and easy use.
- Highly puncture resistant and reliable Polyurethane balloon material for multiple dilatations.
- Enhanced visibility providing higher accuracy in positioning.
- Low profile offering easy advancement through small introducers sheaths / delivery systems.
- Combines perfectly with state-of -the-art stent-grafts.








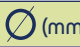
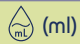


Technical Information






Balloon Material	Polyurethane Material
Balloon Compliance	Compliant
Shaft Diameter	10.5 F
Balloon Marker Distance	40 mm
Recommended max Diameter	46 mm
Balloon Length	46 mm
Rated Burst	50 mm at 100 ml
Shaft Length	65 cm / 100 cm
Recommended Guide Wire	0.035" (0.89 mm)
Recommended Sheath Introducer	12F

!! We recommend the use of an extra stiff guide wire such as the Lokum 'L-Quest' Guide Wire or the Lokum 'Amplatz-Type' Guide Wire.

Sizing Oriented Filling Volume


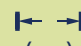
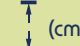

							
 (mm)	10	20	25	30	35	40	46
 (ml)	2	9	15	25	35	45	65

Filling Volume Oriented Sizing

							
 (ml)	5	10	20	30	40	50	65
 (mm)	15	18	27	33	38	42	46

Order Codes

A choice of various diameters and lengths.

Product Code	 (mm)	 (mm)	 (cm)	 (inch)
AL-46-65	46	46	65	0.035"
AL-46-100	46	46	100	0.035"

Office & Warehouse

